



POLICE DEPARTMENT  
**LAKE HAVASU CITY**  
2360 McCULLOCH BOULEVARD NORTH  
LAKE HAVASU CITY, ARIZONA 86403-5947

To Whom It May Concern:

I hereby attest that Josh Sondrol is a Certified Intoxilyzer Quality Assurance Specialist for the Lake Havasu City Police Department.

I certify that the attached documentation of Standard Quality Assurance Procedures and/or Standard Calibration Check Procedure are full, true and correct copy of the original.

*Josh Sondrol*  
Signature

STATE OF ARIZONA        )  
  ) ss  
COUNTY OF MOHAVE    )


Subscribed and sworn to before me this 28 day of MARCH.

2016 by Josh Sondrol

*Carla Duarte*  
NOTARY PUBLIC

3/28/16  
DATE

Commission expires:

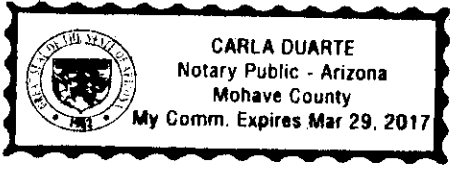

Area Code 928

Administration ..... 855-4884  
Business Office ..... 855-1171  
Investigations ..... 855-5775  
Patrol ..... 855-0515

Fax #'s  
Administration ..... 680-5430  
Business Office ..... 680-5431  
Investigations/Patrol ..... 680-5432

Emergency ..... 9-1-1  
Non-Emergency Dispatch/ ..... 855-4111  
  TDD 855-4114  
Havasu Silent Witness ..... 854-TIPS

E-mail: [police@lhcaz.gov](mailto:police@lhcaz.gov)

EXHIBIT G-2  
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST J SANDER 187 AGENCY LHC PD  
DATE 3/22/16 TIME 2355  
INTOXILYZER SERIAL # 80-003355 LOCATION MAIRJ JAIL

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution .100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"  
(✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.  
(✓) 4. Air blank completed.  
(✓) 5. Calibration check completed. Test results 0.100 AC.  
(✓) 6. Air blank completed.  
(✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE J Sander 187

INTOXILYZER 8000  
Location: MAIN JAIL  
Serial Number: 80-003355  
Core Version: 8105.48  
03/22/2016 23:55:32

-----  
Standard Lot#: 20314100A3  
Last Changed By: J. SONDROL #187  
-----

QAS: J. SONDROL #187  
LAKE HAVASU CITY PD  
-----

Operator: J. SONDROL #187  
LAKE HAVASU CITY PD  
-----

Test	g/210L	Time
Air Blank	0.000	23:56:12
0.100 Cal Check	0.100	23:56:33
Air Blank	0.000	23:57:02

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-4  
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST J SANDROK 187 AGENCY LHC PD  
DATE 3/23/16 TIME 0003  
INTOXILYZER SERIAL # 80-003355 LOCATION MAIN JAIL

(  ) 1. Display Reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.  
(  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.  
(  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(  ) 4. Standard Calibration Check standard 0. 100 AC. Result 0. 100 AC.

Instrument is operating properly and accurately. YES  NO

COMMENTS

SIGNATURE J Sandrok 187

INTOXILYZER 8000  
Location: MAIN JAIL  
Serial Number: 80-003355  
Core Version: 8105.48  
03/23/2016

00:03:50

---

Standard Lot#: 20314100A3  
Last Changed By: J. SONDROL #187

---

QAS: J. SONDROL #187  
LAKE HAVASU CITY PD

---

Operator: J. SONDROL #187  
LAKE HAVASU CITY PD

---

Subject: SQAP,TEST  
DOB: 01/01/1950  
Sex: M Weight: 123

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	00:04:39
Diagnostic Test	Pass	00:05:14
Air Blank	0.000	00:05:41
0.100 Cal Check	0.100	00:06:02
Air Blank	0.000	00:06:31
Subject Test	0.000	00:06:49
Air Blank	ABT*	00:06:56
Air Blank	IPS**	00:07:03

\*Sequence Aborted

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: MAIN JAIL  
Serial Number: 80-003355  
Core Version: 8105.48  
03/23/2016

00:07:45

---

Standard Lot#: 20314100A3  
Last Changed By: J. SONDROL #187

---

QAS: J. SONDROL #187  
LAKE HAVASU CITY PD

---

Operator: J. SONDROL #187  
LAKE HAVASU CITY PD

---

Subject: SQAP,TEST  
DOB: 01/01/1950  
Sex: M Weight: 123

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	00:08:30
Diagnostic Test	Pass	00:09:05
Air Blank	0.000	00:09:33
0.100 Cal Check	0.100	00:09:53
Air Blank	0.000	00:10:23
Subject Test	0.203*	00:13:30
Air Blank	0.000	00:14:03
Air Blank	0.000	00:14:31
Subject Test	XXX**	00:14:55
Air Blank	RFI***	00:15:03

\*Deficient Sample  
\*\*Invalid Sample  
    Begin new deprivation period  
\*\*\*RFI Detect

**Not a Successfully Completed Test Sequence**

EXHIBIT G-2  
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST J SANDRA 187 AGENCY LNC PD  
DATE 3/23/16 TIME 0016  
INTOXILYZER SERIAL # 80-003355 LOCATION MAIN JAIL

- (  ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution .100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (  ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (  ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (  ) 4. Air blank completed
- (  ) 5. Calibration check completed. Test results 0. 100 AC.
- (  ) 6. Air blank completed.
- (  ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE C. Sandra 187

INTOXILYZER 8000  
Location: MAIN JAIL  
Serial Number: 80-003355  
Core Version: 8105.48  
03/23/2016

00:16:13

-----  
Standard Lot#: 20314100A3  
Last Changed By: J. SONDROL #187  
-----

QAS: J. SONDROL #187  
LAKE HAVASU CITY PD  
-----

Operator: J. SONDROL #187  
LAKE HAVASU CITY PD  
-----

Test	g/210L	Time
Air Blank	0.000	00:16:42
0.100 Cal Check	0.100	00:17:03
Air Blank	0.000	00:17:33

Preventative Maintenance  
Performed -- Timer Reset